

SURFLO™ HYBRIA™ Closed I.V. Catheter
Instruction for use

TERUMO MEDICAL CARE SOLUTIONS

Needle & Catheter Insertion

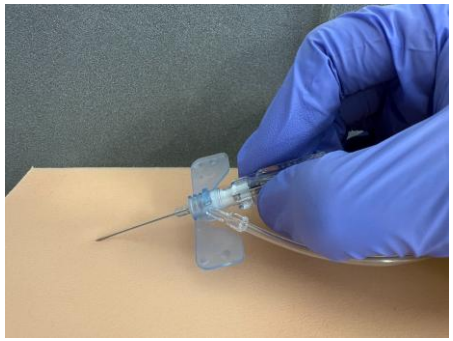


Hold the needle hub with needle bevel facing upwards. Then insert. Confirm that blood returns along the catheter and catheter tip is inside the vessel.

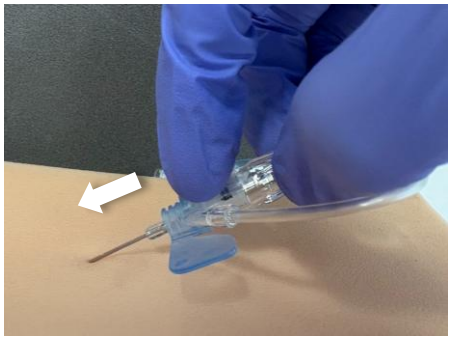


Stabilize the needle hub with finger, advance the catheter to desired depth.

SurFlash™


Visual confirmation of catheter placement occurs due to a unique flashback feature.

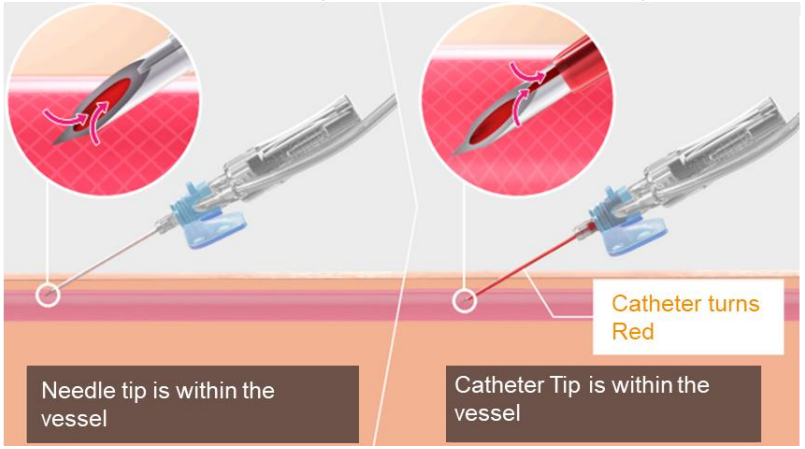

Needle withdrawal

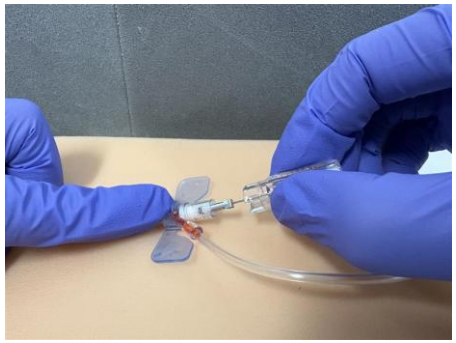

Stabilize the catheter hub and tab with finger, withdraw needle in a straight manner until safety shield is fully activated.

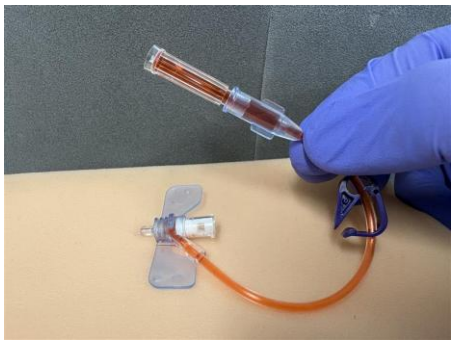


Then blood return should be observed in the extension tube.

Catheter securement and I.V. line connection



Apply a film dressing to protect the puncture site.

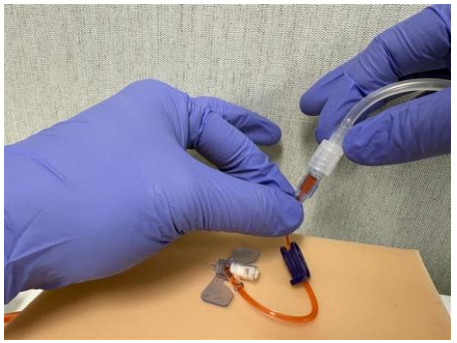


Close the clamp and remove vent fitting(s). Connect the extra needle-free connector provided in the primary package or other infusion set.

SURFLO™ HYBRIA™ Closed I.V. Catheter
Instruction for use


The insertion by holding the wings

Needle & Catheter Insertion

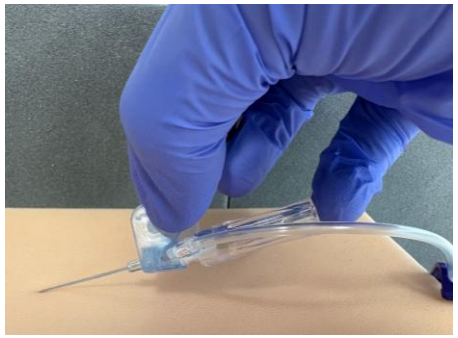

Fold the wings and hold with needle bevel facing upwards. Then insert. Confirm that blood returns along the catheter and catheter tip is inside the vessel.


Note : Put ring finger or little finger behind needle hub to prevent needle from retracting during insertion.

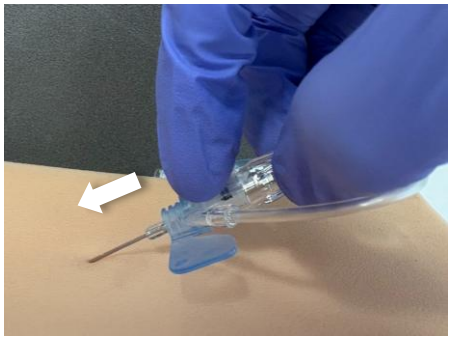


Stabilize the needle hub with finger, advance the catheter to desired depth. Stabilize the catheter hub and tab with finger, withdraw needle in a straight manner until safety shield is fully activated.

PRECAUTIONs

During needle withdrawal

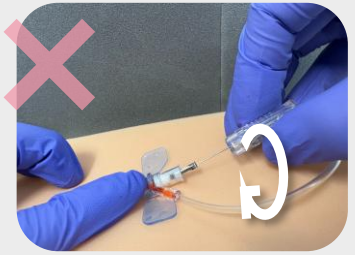

Do not withdraw the needle while rotating

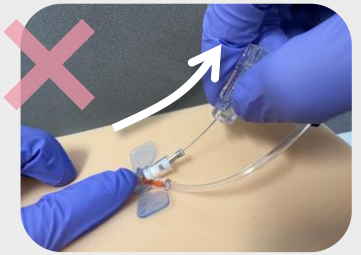

Do not withdraw the needle at an angle

〔To prevent safety shield deactivation〕


During Line connection



Make sure to close Clamp before vent fitting removal

〔To prevent blood exposure〕

During connection with needle-free connector

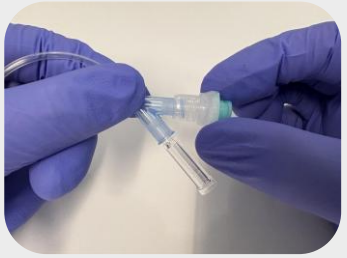

Make sure that the connection of injection port is not loose

RX ONLY. This advertisement is directed to physicians only, and not to consumers.
Refer to product labels and packaging insert for complete warnings, precautions, potential complications, and instructions for use.

TERUMO MEDICAL CARE SOLUTIONS

@2025 TERUMO MEDICAL CORPORATION
265 Davidson Avenue Somerset NJ 08873 U.S.A